




**WORLD
PRECISION
INSTRUMENTS**
Instrumenting scientific ideas

SPL Syringe Pumps

Superior reliability and flow performance

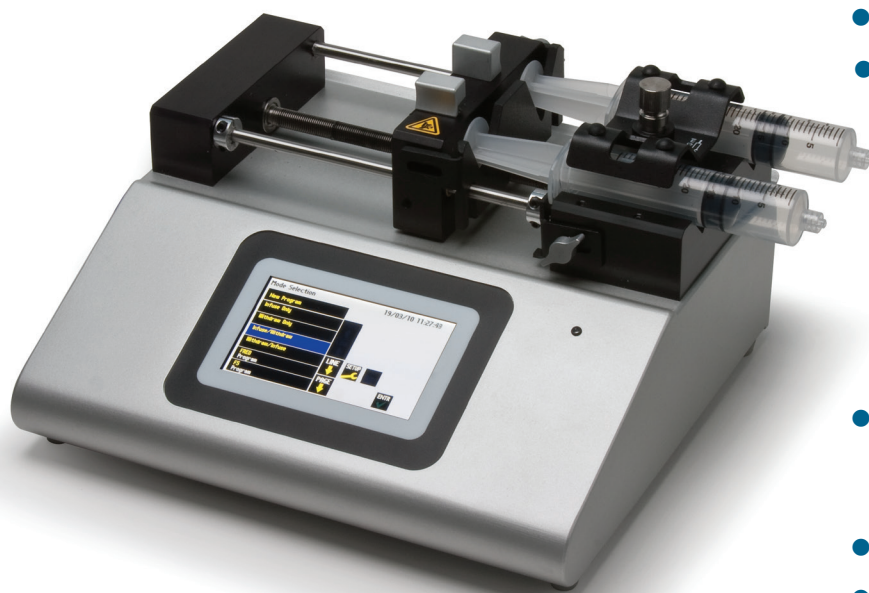

- *Easy to configure*
- *Easy to use*
- *Vertical mounting*
- *Quick release pusher block*
- *Wide clamp for larger syringes*
- *Multi-size syringes*
- *Anti-siphon bracket*
- *RS232 computer connection, USB serial input, RS485 for networking pumps*
- *Wide variety of applications*
- *Worldwide regulatory compliance*

This new SPL series of pumps show all the pertinent, real-time information on a single touch screen display. Information displayed includes flow rate, grading syringe volume, total remaining time and total volume delivered. Recognizable, international icons make it easy to use.

To conserve table space, mount the SPL pump vertically. The display automatically reorients itself for easy vertical operation.

The welded, unibody, steel chassis of the SPL pumps makes them quieter to operate. It also offers less vibration, no deformation and excellent EMI/RFI shielding.

The one-touch operation makes the

SPL pump easy to configure, and its 15-character labeling system makes saved programs quick to identify. You can easily transfer recipes between pumps for consistent configuration without reprogramming. Pumps can be connected for gradients, increased flow capacity and simultaneous multi-tests.

This family of pumps includes:

- Infuse only (WPI #**SPLG200**)
- Infuse/withdraw (WPI #**SPLG210**)
- Infuse/withdraw programmable (WPI #**SPLG212**)
- Push-pull (WPI #**SPLG270**)
- Push-pull programmable (WPI #**SPLG272**)

SPL Syringe Pumps

Superior reliability and flow performance



When mounted vertically, the display screen of the SPL series pumps automatically reorients for ease of use.

SPECIFICATIONS

Accuracy	±0.35%
Reproducibility	±0.05%
Syringes (min./max.)	0.5µL/140mL
Flow rate minimum (0.5µL syringe)	5pL/min.
Flow rate maximum (140mL syringe)	220.97mL/min.
Display	4.3" WQVGA TFT color display with touch pad
Non-volatile memory	Stores all settings

Connectors:

RS-232	9-pin D-sub connector
RS-485	IEEE-1394, 6 pos
USB	type B
I/O and TTL	15-pin D-sub connector
Linear force (maximum)	34kg (75 lb.) @ 100% force selection
Drive motor	1.8° stepper motor
Motor drive control	Microprocessor with 1/16 microstepping
Number of microsteps per one rev. of lead screw	6,400
Step rate minimum	27.5 sec./µstep
Step rate maximum	26µsec./µstep
Pusher travel rate minimum	0.36µm/min.
Pusher travel rate maximum	190.80mm/min.
Power	100-240VAC; 50/60Hz; 50W; 0.5A fuse
Dimensions	8.89 x 25.4 x 27.94cm (3.5 x 10 x 11")
Weight	4.9kg (10.75 lb.)

Atmospheric specifications:

Operating temperature	4–40°C (40–104°F)
Storage temperature	–10–70°C (14–158°F)
Humidity	20–80%RH, non-condensing
Mode of operation	Continuous
Classification	Class I
Pollution degree	1
Installation Category	II
Regulatory Certifications	CE, UL, CSA, CB Scheme, EU RoHS



WORLD PRECISION INSTRUMENTS

USA: International Trade Center, 175 Sarasota Center Boulevard, Sarasota FL 34240-9258 USA

Tel: 941-371-1003 • **Fax:** 941-377-5428 • **E-mail:** wpi@wpiinc.com • **Internet:** www.wpiinc.com

UK: 1 Hunting Gate, Hitchin, Hertfordshire SG4 0TJ England • Tel: 44 (0)1462 424700 • E-mail: wpiuk@wpi-europe.com

Germany: Zossener Str. 55, 10961 Berlin, Germany • Tel: 030-6188845 • Fax: 030-6188670 • E-mail: wpide@wpi-europe.com

China & Hong Kong: Rm 29a, No8 Donfang Rd., Pudong District, Shanghai 200120 PRC • Tel: +86 688 85517 • E-mail: ChinaSales@china.wpiinc.com

Brazil: Conselheiro Nabias, 756 sala2611, Santos-Sao Paulo 11045-002 Brazil • E-mail: info@brazil.wpiinc.com